


FALL

COLOUR GUIDE

LARYVOX® MY EXTRA HME

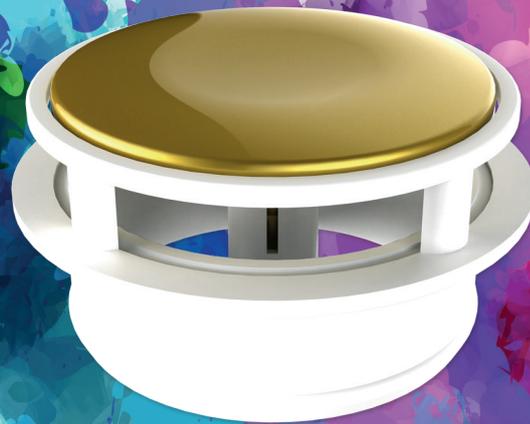

MAKE YOUR LIFE MORE COLOURFUL -

YOUR PERSONAL HME FILTER

FOREWORD

MAKE YOUR LIFE MORE COLOURFUL WITH THE NEW LARYVOX® MY EXTRA HME

You are laryngectomised and regularly wear HME filters? You are brave and do not hide your tracheostoma? The new LARYVOX® MY EXTRA HME filter is definitely the right choice for all self-confident laryngectomised patients! Test our colourful MY EXTRA HME.

The product and filter properties are identical to those of the classic LARYVOX® EXTRA HME. The only difference is the colour.

FILTERING - WARMING - HUMIDIFYING

The regular use of HME filters is recommended to compensate for the loss of nasal respiration. This filter is used for heat and moisture exchange and filters the breathing air to prevent irritation of the airways due to dry air, dust or harmful substances contained in the breathing air. After surgery, tracheotomised/laryngectomised patients frequently suffer from increased production of mucous and coughing. These symptoms can impair breathing and voice rehabilitation. Regular use of HME filter cassettes alleviates the symptoms through reduction of viscous secretions in the airways and enhances the quality of life.

The filter cassette is intended for single use only and must be replaced at the latest after 24 hours. LARYVOX® MY EXTRA HME can be worn both during the day and at night.

EASY SPEAKING OPTION FOR PATIENTS WITH SHUNT VALVES

The LARYVOX® MY EXTRA HME is equipped with an easy-to-operate speaking option for patients with a shunt valve.

The rounded recess in the lid of LARYVOX® MY EXTRA HME additionally facilitates deliberate closure of the tracheostoma and renders the optimal finger position easy to feel. By exerting slight finger pressure on the HME lid, the patient is able to accomplish a very precise closure of the tracheostoma, for instance for speaking. The device can be operated even through clothing covering it.

FILTERS TO MEET YOUR NEEDS

By virtue of four different versions, the LARYVOX® MY EXTRA HME can be selected to suit individual activities and needs:

- Normal** ▶ normal breathing resistance for situations requiring normal physical activity
- Medium** ▶ medium breathing resistance for physical activities
- HighFlow** ▶ reduced breathing resistance for situations requiring strong physical exertion
- Sport** ▶ very low respiratory resistance specifically for sporting activities



CONFIGURATION

EASE OF USE

Lateral breathing apertures enable largely unhindered air flow even if the cassette is covered with clothing or additional tracheostoma protection material, thus enabling easy breathing. The HME filter cassette can be attached to all conventional base plates (e.g. LARYVOX® Tape) or silicone tracheostomy tubes for laryngectomised patients (e.g. LARYNGOTEC®) with a 22 mm adapter (KOMBI). A finger ledge on the housing enables easy handling and quick removal of the LARYVOX® MY EXTRA HME filter cassette from the holding ring, thereby facilitating expectoration and tracheostoma care.

NO UNPLEASANT BACKGROUND NOISES

The HME lid and housing are optimally matched to each other through a special closure mechanism. The HME is thus sealed particularly tightly, preventing unintentional air loss from the side of the HME. Unpleasant side noises are reduced.

YOUR OWN DESIGN

You can now create your own personal MY EXTRA HME. To do this, first select your filter: Normal, Medium, High-Flow or Sport. You can then determine the colour of the housing and the speech button (lid) and mix and match the colours.

You can choose from 23 different colours and over 500 different colour combinations. The colours of the cover and housing can be combined in any combination. This gives you the option to choose your HME daily to go with your clothes.

HAVE WE AROUSED YOUR INTEREST?

Visit us at www.fahl.de and discover the colour variety of the LARYVOX® MY EXTRA HME filter. The service employee responsible for you in the office and the field force representative of our company will be pleased to help you in designing your own HME filters!

Simply call us!

HOW TO DESIGN YOUR PERSONAL HME

1. Select the LARYVOX® MY EXTRA HME filter:

- Normal (REF 49860-)
- Medium (REF 49862-)
- HighFlow (REF 49861-)
- Sport (REF 49863-)

2. Choose the housing colour:


- e.g. Silver (colour code 99)

3. Choose colour for lid/pressure piece:

- e.g. Turquoise (colour code 65)

Example:

- ① Filter: HighFlow REF 49861-
- ② Colour housing: Silver -99
- ③ Colour lid: Turquoise -65
▶ REF 49861-9965

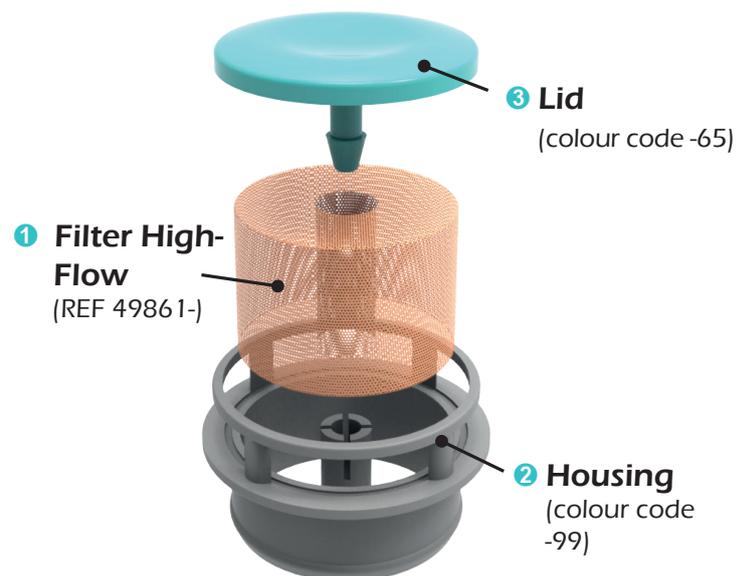

LARYVOX® MY EXTRA HME
REF 49861-9965

COLOUR SHADES



NATURAL SHADES



EXAMPLE COLOUR COMBINATIONS

LARYVOX® MY EXTRA HME



Classic/Yellow



Classic/Gold



Classic/Brown



Classic/Orange



Classic/Red



Classic/Purple



Classic/Blue



Classic/Turquoise



Classic/Green



Classic/Black



Classic/Silver



Classic/White



Yellow/Yellow



Yellow/Red



Yellow/Dark blue



Yellow/Green



Yellow/Black



Yellow/White



Gold/Gold



Gold/Silver



Gold/White



Brown/Brown



Brown/White



Orange/Orange



Orange/White



Red/Yellow



Red/Red



Red/Black



Red/Silver



Red/White



Purple/Purple



Purple/Black



Purple/Silver



Purple/White



Dark blue/Yellow



Dark blue/Dark blue



Dark blue/White



Blue/Blue



Blue/Green



Blue/Black

EXAMPLE COLOUR COMBINATIONS

LARYVOX® MY EXTRA HME


| | | | | |
|----------------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------------------------------|
|  <p>5599</p> <p>Blue/Silver</p> |  <p>5590</p> <p>Blue/White</p> |  <p>6565</p> <p>Turquoise/Turquoise</p> |  <p>6599</p> <p>Turquoise/Silver</p> |  <p>6590</p> <p>Turquoise/White</p> |
|  <p>6010</p> <p>Green/Yellow</p> |  <p>6060</p> <p>Green/Green</p> |  <p>6090</p> <p>Green/White</p> |  <p>9510</p> <p>Black/Yellow</p> |  <p>9598</p> <p>Black/Gold</p> |
|  <p>9530</p> <p>Black/Red</p> |  <p>9545</p> <p>Black/Purple</p> |  <p>9555</p> <p>Black/Blue</p> |  <p>9565</p> <p>Black/Turquoise</p> |  <p>9560</p> <p>Black/Green</p> |
|  <p>9595</p> <p>Black/Black</p> |  <p>9599</p> <p>Black/Silver</p> |  <p>9590</p> <p>Black/White</p> |  <p>9910</p> <p>Silver/Yellow</p> |  <p>9998</p> <p>Silver/Gold</p> |
|  <p>9945</p> <p>Silver/Purple</p> |  <p>9955</p> <p>Silver/Blue</p> |  <p>9965</p> <p>Silver/Turquoise</p> |  <p>9995</p> <p>Silver/Black</p> |  <p>9999</p> <p>Silver/Silver</p> |
|  <p>9990</p> <p>Silver/White</p> |  <p>9010</p> <p>White/Yellow</p> |  <p>9098</p> <p>White/Gold</p> |  <p>9080</p> <p>White/Brown</p> |  <p>9020</p> <p>White/Orange</p> |
|  <p>9030</p> <p>White/Red</p> |  <p>9045</p> <p>White/Purple</p> |  <p>9050</p> <p>White/Dark blue</p> |  <p>9055</p> <p>White/Blue</p> |  <p>9065</p> <p>White/Turquoise</p> |
|  <p>9060</p> <p>White/Green</p> |  <p>9095</p> <p>White/Black</p> |  <p>9099</p> <p>White/Silver</p> |  <p>9090</p> <p>White/White</p> | |



ORDER INFORMATION:

- | | | |
|----------------------------------|-----------|----------------|
| ▶ LARYVOX® MY EXTRA HME | Pack of 5 | REF 49860-xxxx |
| ▶ LARYVOX® MY EXTRA HME MEDIUM | Pack of 5 | REF 49862-xxxx |
| ▶ LARYVOX® MY EXTRA HME HIGHFLOW | Pack of 5 | REF 49861-xxxx |
| ▶ LARYVOX® MY EXTRA HME SPORT | Pack of 5 | REF 49863-xxxx |

When ordering, please quote colour code after the article number (REF)!

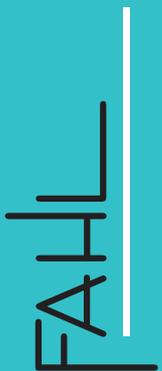

Andreas Fahl
Medizintechnik-Vertrieb GmbH
August-Horch-Str. 4a
51149 Cologne - Germany
Phone +49 (0) 22 03/29 80-0
Fax +49 (0) 22 03/29 80-100
mail info@fahl.de
www.fahl.de

Local branch Berlin
Am Borsigturm 62
13507 Berlin - Germany
Phone +49 (0) 30/43 73 139-0
Fax +49 (0) 30/43 73 139-35

Local branch Munich
Aurbacher Straße 2
81541 Munich - Germany
Phone +49 (0) 89/58 90 90 9-0
Fax +49 (0) 89/58 90 90 9-20

Overseas branch Austria
Triester Straße 14 / Bürohaus 2
2351 Wiener Neudorf - Austria
Phone +43 (0) 22 36/89 28 38-0
Fax +43 (0) 22 36/89 28 38-14

Overseas branch Switzerland
Cherstrasse 4
8152 Glattbrugg - Switzerland
Phone +41 (0) 44829360
Fax +41 (0) 48293061