
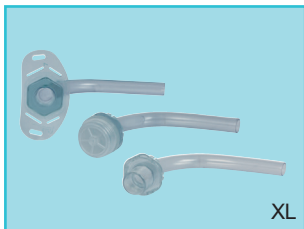


DURAVENT® UNI-KOMBI LINGO-PHON

1 ICU + 1 ICFK +
1 COMBIPHON®

- ▶ Tracheostomy tube for long-term use, can be used up to 3-6 months
- ▶ With two inner cannulas (IC): One fenestrated with 22 mm adapter (ICFK / KOMBI) to connect the enclosed COMBIPHON® speaking valve (closed position valve with silicone membrane), one closed with 15 mm connector (ICU / UNI) to connect an HME
- ▶ Easy-Lock has been designed for fixing the inner cannula and compatible accessories securely - no twisting
- ▶ Outer cannula (AK) with sieve-fenestration
- ▶ Two length variants: standard and XL
- ▶ Soft, pliable, lightweight
- ▶ Transparent, inconspicuous design
- ▶ Custom-made variants possible (length, sieve-fenestration, bending)
- ▶ Including broad adjustable OPTIFLAUSCH® tube holder and a product ID card

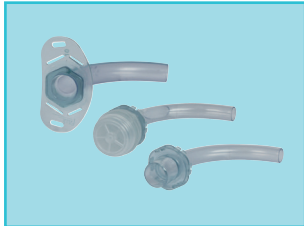

Order information:

DURAVENT® UNI-KOMBI LINGO-PHON, SIZE 5-13



REF 21072

1 ICU + 1 ICFK + 1 COMBIPHON®

DURAVENT® UNI-KOMBI LINGO-PHON XL, SIZE 5-13



REF 21272

1 ICU + 1 ICFK + 1 COMBIPHON®

When ordering, please state the size after the item number!

Size	O.D. neck flange mm	O.D. Tip mm	I.D. mm O-cannula	I.D. mm I-cannula	Length mm standard	Length mm extra-long (XL)
5	8.4	6.6	5.4	4.3	55	95
6	9.4	7.8	6.2	4.9	60	95
7	10.2	8.8	7.0	5.5	65	100
8	11.2	9.8	8.0	6.5	70	105
9	12.4	11.2	9.0	7.5	75	105
10	13.4	12.2	10.0	8.5	85	110
11	14.4	13.2	11.0	9.5	87	115
12	15.4	14.2	12.0	10.2	90	120
13	16.6	15.2	13.0	11.2	90	125

O.D. neck flange = Outer diameter behind the neck flange · O.D. Tip = Outer diameter of the outer cannula at the tip of the cannula · I.D. = Inner diameter at the tip of the cannula · O-cannula = Outer cannula · I-cannula = Inner cannula · Length = Length over the outer curve



Keep away from sunlight



Single patient use



Consult instructions for use

CAUTION:
United States Federal law restricts this device to sale, distribution and use by or on order of a physician or a licensed practitioner (in the US only).



CE marking identification number of the notified body

Product made in Germany