
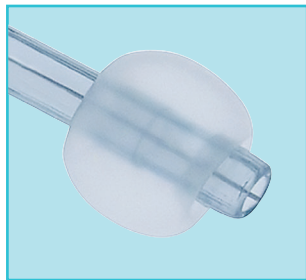


# TRACHEOTEC® VARIO CUFF-SUCTION

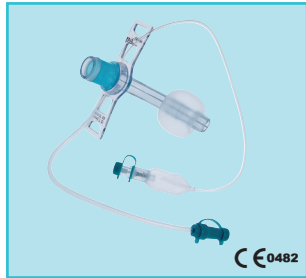

- ▶ Tracheostomy tube with cuff and integrated suction line to remove secretions above the CUFF
- ▶ Suction line integrated into the inner wall of the tube to avoid any additional irritation of the trachea or the areas around the tracheostoma
- ▶ Soft, flexible neck-flange with key tube identification data clearly marked
- ▶ With rotating 15 mm connector (VARIO) to buffer torque and tension caused by tractive forces from ventilation equipment
- ▶ With X-ray contrast stripe
- ▶ Set includes tube holder, obturator and product ID card

## Order information:

**TRACHEOTEC® VARIO CUFF-SUCTION SIZE 5-10**

When ordering, please state the size after the item number!



**REF 18930**

Size	O.D. neck flange (mm)	O.D. Tip (mm)	I.D. (mm)	Length midline (mm)	Length (mm)	REF
5	7.3	7.3	5.0	55.7	61.1	<b>18930-05</b>
6	8.7	8.7	6.0	60.7	67.1	<b>18930-06</b>
7	10.0	10.0	7.0	69.6	77.0	<b>18930-07</b>
8	11.3	11.3	8.0	79.5	87.9	<b>18930-08</b>
9	12.7	12.7	9.0	85.5	94.9	<b>18930-09</b>
10	14.0	14.0	10.0	95.0	105.3	<b>18930-10</b>

O-tube = outer cannula · O.D. neck flange = outer diameter behind the neck flange · O.D. Tip = outer diameter at the tube tip · I.D. = inner diameter at cannula tip · I-tube = inner cannula · Length = length over the outer curve · Length at mid radius = length of centre line from neck flange to tip