


# BLOM-SINGER® DUAL VALVE™ INDWELLING / DUAL VALVE™ INDWELLING LF (LARGE FLANGE) VOICE PROSTHESIS

- ▶ The Blom-Singer® DUAL VALVE™ Indwelling prosthesis is equipped with two independent consecutive low-pressure valves made of silver oxide-containing silicone
- ▶ Should one of the two valves fail due to candida infestation, then the second valve ensures that proper functioning is maintained
- ▶ This Blom-Singer® DUAL VALVE™ Indwelling LF (Large Flange) voice prosthesis offers a functional extension compared to the Blom-Singer® DUAL VALVE™ Indwelling
- ▶ It is additionally equipped with two enlarged flanges (oesophageal and tracheal flange) for treatment of periprosthetic leakages as appropriate for the indication
- ▶ Both voice prosthesis are available as 20 Fr. version in sizes 6 - 14
- ▶ With the help of the LARYVOX® INSERTER the voice prosthesis can be placed anterograde quickly and easily:  
 LARYVOX® INSERTER (REF 25900) - for different types of voice prosthesis  
 LARYVOX® INSERTER EXPERT (REF 25910) - for 16 Fr. voice prosthesis (or DUAL VALVE™ 20 Fr.)  
  
 Application with gel capsule technique:  
 LARYVOX® INSERTER SPECIAL - (REF 25920) - for difficult tracheostoma

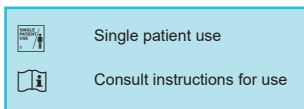

### Order information:

Blom-Singer® DUAL VALVE™ Indwelling	20 Fr.	SIZE 6-14	REF DV20XX-NS
Blom-Singer® DUAL VALVE™ Indwelling LF (Large Flange)	20 Fr.	SIZE 6-14	REF DV20XX-LF
LARYVOX® INSERTER			REF 25900
LARYVOX® INSERTER EXPERT			REF 25910
LARYVOX® INSERTER SPECIAL			REF 25920