

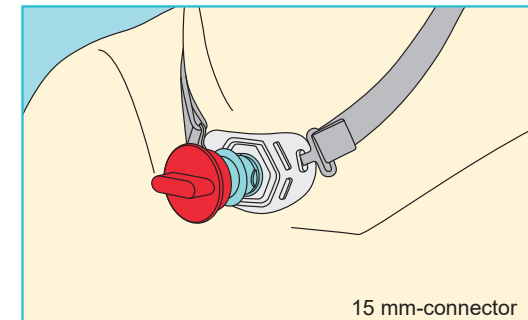


FAHL® MULTI PLUG

Decannulation plug

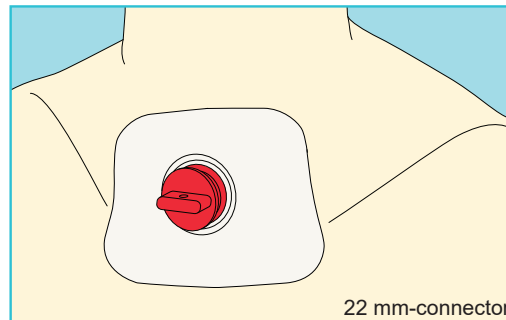

- ▶ The FAHL® MULTI PLUG serves to seal a tracheostomy tube/base plate and enables safe decannulation
- ▶ After a tracheotomy, the patient must be slowly re-accustomed to breathing normally via the physiological airway with a temporary tracheostoma closure
- ▶ During application the patient is ventilated solely via the nasal and pharyngeal cavity
- ▶ Decannulation can be performed when the primary cause of the tracheostoma is no longer present, i.e. in tracheotomised patients with partially or completely preserved larynx
- ▶ The FAHL® MULTI PLUG has a bar for easier handling during insertion/removal
- ▶ An opening integrated in the bar allows the attachment of a strap to protect against loss, e.g. in case of a coughing fit
- ▶ The FAHL® MULTI PLUG is compatible with all common 22 mm, e.g. LARYVOX® Tape base plate and 15 mm-mounting systems, e.g. SPIRAFLEX® UNI
- ▶ The FAHL® MULTI PLUG cannot be used in conjunction with TRACHEOTEC®, TRACHEOTEC® PRO/ SILC and TRACHEOSILC®

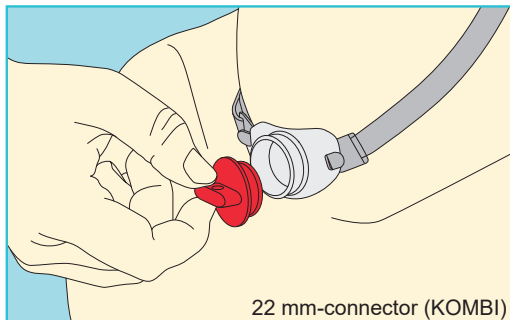

Order information:
FAHL® MULTI PLUG

REF 49812