
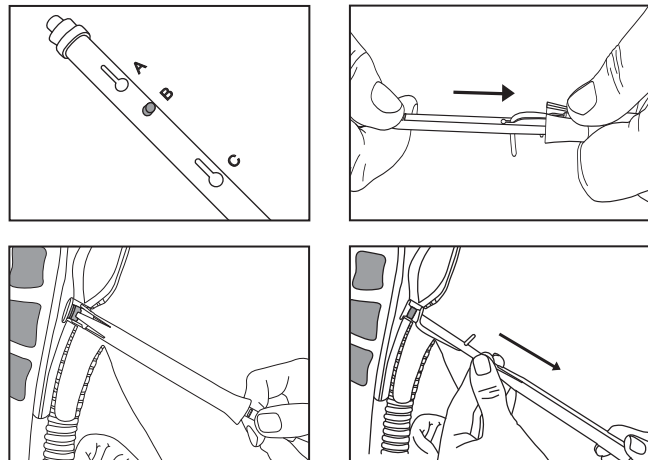


# LARYVOX® INSERTER

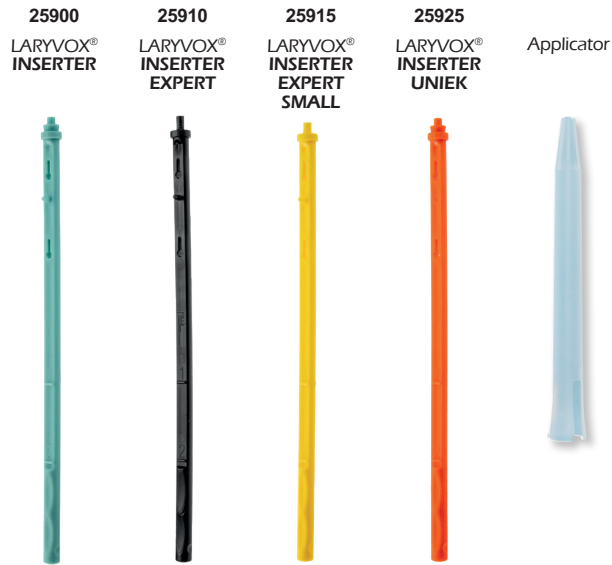

- ▶ The **LARYVOX® INSERTER** consists of a placement instrument as well as an adapter and facilitates the perfect anterograde placement or replacement of a shunt valve
  - **REF 25900**: for voice prostheses with 5-14 mm shaft length
  - **REF 25910**: for a 16 Fr. voice prosthesis and/or a Blom-Singer® DUAL VALVE shunt valve (20 Fr.)
  - **REF 25915**: for voice prostheses with a shaft length of 4 mm
  - **REF 25925**: for **LARYVOX® UNIEK** voice prosthesis/shunt valve
- ▶ The placement instrument has two (**REF 25925**) resp. three different fixation options to allow easy and safe placement of various types of shunt valves.
- ▶ The applicator of the **LARYVOX® INSERTER** enables effortless insertion of the shunt valve and is readily prepared in a few steps.
- ▶ Additionally, the rounded applicator tip promotes gentle placement.

## SCOPE OF DELIVERY

- ▶ Applicator
- ▶ Instruction for use

**LARYVOX® INSERTER** REF 25900  
1 x placement instrument (turquoise)

or  
**LARYVOX® INSERTER EXPERT** REF 25910  
1 x placement instrument (black)

or  
**LARYVOX® INSERTER EXPERT SMALL** REF 25915  
1 x placement instrument (yellow)

The voice prosthesis is not included in delivery and must be ordered separately. Please heed the instructions for use of the respective voice prosthesis manufacturer!

The sets of the voice prosthesis **LARYVOX® UNIEK** (**REF 25010** und **25020**) contain the:

**LARYVOX® INSERTER UNIEK** REF 25925  
1 x placement instrument (orange)


Product	Unit	REF
<b>LARYVOX® INSERTER</b>	1 piece	<b>25900</b>
<b>LARYVOX® INSERTER EXPERT</b>	1 piece	<b>25910</b>
<b>LARYVOX® INSERTER EXPERT SMALL</b>	1 piece	<b>25915</b>
<b>LARYVOX® INSERTER UNIEK</b>	1 piece	<b>25925</b>



Andreas Fahl Medizintechnik-Vertrieb GmbH is the exclusive distributor for Blom-Singer® products in selected European countries. Blom-Singer® is a registered trademark of Freudenberg Medical, LLC.



# LARYVOX® INSERTER

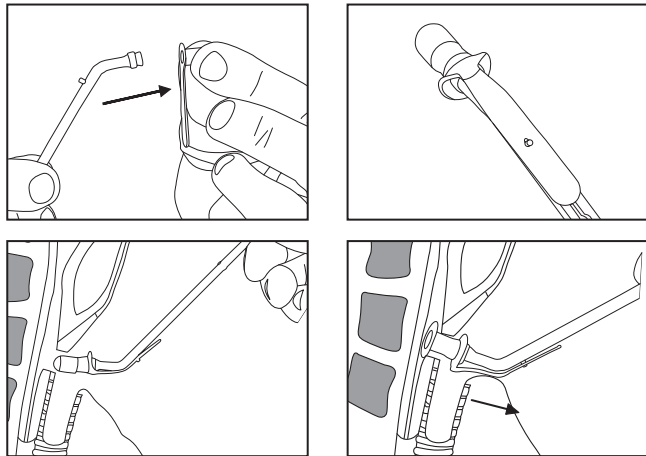

- ▶ **LARYVOX® INSERTER SPECIAL** is used for optimal antegrade placement or replacement of a shunt valve by application of the gel capsule technique in cases of:
  - deep-positioned shunt or shunt with oblique course
  - constricted shunt channel
  - particularly small tracheostoma for optimising the viewing conditions
- ▶ The placement instrument has one safe fixation pin to allow easy and safe placement of various shunt valve types with a diameter of 20 French from 5 mm in length
- ▶ Effortless insertion of the shunt valve is enabled by application of the gel capsule technique
- ▶ The shunt valve is introduced atraumatically into the shunt channel with the aid of the gel capsule
- ▶ A drop of water-soluble local anaesthetic can be applied directly prior to placement in order to accelerate the dissolution of the gel capsule in the oesophagus
- ▶ The matching gel capsules are available separately in sizes 20 Fr. and 22 Fr.

## SCOPE OF DELIVERY

- ▶ Instruction for use

**LARYVOX® INSERTER SPECIAL** REF 25920  
1 x placement instrument (blue)

Optional accessories:  
**BLOM-SINGER® REPLACEMENT GEL CAPS** REF BE3X9X  
Replacement Gel Caps

The shunt valve is not included in delivery and must be ordered separately. Please heed the instructions for use of the respective shunt valve manufacturer!

Product	Unit	REF
LARYVOX® INSERTER SPECIAL	1 piece	25920













Andreas Fahl Medizintechnik-Vertrieb GmbH is the exclusive distributor for Blom-Singer® products in selected European countries. Blom-Singer® is a registered trademark of Freudenberg Medical, LLC.










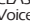
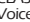
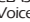
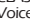
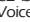

















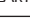


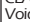














# LARYVOX® INSERTER
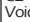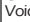
## ASSIGNMENT

25900	25910	25915	25920	25925
LARYVOX® INSERTER	LARYVOX® INSERTER EXPERT	LARYVOX® INSERTER EXPERT SMALL	LARYVOX® INSERTER SPECIAL	LARYVOX® INSERTER UNIEK
Color				
				
				

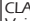Take into consideration beforehand the length (xx) of the voice prosthesis to be inserted.

REF	Name	Lenght	INSERTER
20 FRENCH			
AD20xx-NS	ADVANTAGE Indwelling 20 Fr. Voice Prosthesis	06-14 mm	
AD2004-NS	ADVANTAGE Indwelling 20 Fr. Voice Prosthesis	04 mm	
AS20xx-NS	ADVANTAGE Indwelling 20 Fr. Voice Prosthesis	06-14 mm	
AS2004-NS	ADVANTAGE Indwelling 20 Fr. Voice Prosthesis	04 mm	
DV20xx-LF	DUAL VALVE Indwelling LF 20 Fr. Voice Prosthesis	06-14 mm	
DV20xx-NS	DUAL VALVE Indwelling 20 Fr. Voice Prosthesis	06-14 mm	
IN2004-NS	CLASSIC Indwelling 4 mm 20 Fr. Voice Prosthesis	04 mm	
IN20xx-NS	CLASSIC Indwelling 20 Fr. Voice Prosthesis	06-20 mm	
IN2004-LEF	CLASSIC Indwelling LEF 4 mm 20 Fr. Voice Prosthesis	04 mm	
IN20xx-LEF	CLASSIC Indwelling LEF 20 Fr. Voice Prosthesis	06-14 mm	
IN20xx-IR	CLASSIC Indwelling IR 20 Fr. Voice Prosthesis	06-12 mm	
IN2004-LEIR	CLASSIC Indwelling LEIR 20 Fr. 4 mm Voice Prosthesis	04 mm	
IN20xx-LEIR	CLASSIC Indwelling LEIR 20 Fr. Voice Prosthesis	06-08 mm	
IN20xx-LESL	CLASSIC Indwelling LESL 20 Fr. Voice Prosthesis	05 und 07 mm	
IN2004-LETO	TEP Occluder Indw. LETO 20 Fr. 4 mm Voice Prosthesis	04 mm	
IN20xx-LETO	TEP Occluder Indw. LETO 20 Fr. Voice Prosthesis	06-14 mm	
IN2004-LF	CLASSIC Indwelling LF 20 Fr. 4 mm Voice Prosthesis	04 mm	
IN20xx-LF	CLASSIC Indwelling LF 20 Fr. Voice Prosthesis	06-14 mm	
IN20xx-SL	CLASSIC Indwelling SL 20 Fr. Voice Prosthesis	05,07 und 09 mm	
IN20xx-SO	CLASSIC Indwelling STERIL 20 Fr. Voice Prosthesis	08-14 mm	
25020-xx	LARYVOX® Uniek U.L.R.	05-13 mm	
25010-xx	LARYVOX® Uniek U.L.R.	05-13 mm	

REF	Name	Lenght	INSERTER
16 FRENCH			
AS1604-NS	ADVANTAGE Indwelling 16 Fr. Voice Prosthesis	04 mm	
AS16xx-NS	ADVANTAGE Indwelling 16 Fr. Voice Prosthesis	06-14 mm	
IN1604-NS	CLASSIC Indwelling 4 mm 16 Fr. Voice Prosthesis	04 mm	
IN16xx-NS	CLASSIC Indwelling 6 mm 16 Fr. Voice Prosthesis	06-20 mm	
IN16xx-SO	CLASSIC Indwelling STERIL 10 mm 16 Fr. Voice Prosthesis	08-14 mm	
IN16xx-IR	CLASSIC Indwelling IR 16 Fr. Voice Prosthesis	06-12 mm	
IN1604-LEF	CLASSIC Indwelling LEF 4 mm 16 Fr. Voice Prosthesis	04 mm	
IN16xx-LEF	CLASSIC Indwelling LEF 16 Fr. Voice Prosthesis	06-14 mm	
IN16xx-LESL	CLASSIC Indwelling LESL 16 Fr. Voice Prosthesis	05 und 07 mm	
IN1604-LETO	CLASSIC Indwelling LETO 4 mm 16 Fr. Voice Prosthesis	04 mm	
IN16xx-LETO	CLASSIC Indwelling LETO 16 Fr. Voice Prosthesis	06-14 mm	
IN1604-LF	CLASSIC Indwelling LF 4 mm 16 Fr. Voice Prosthesis	04 mm	
IN16xx-LF	CLASSIC Indwelling LF 16 Fr. Voice Prosthesis	06-14 mm	
IN16xx-SL	CLASSIC Indwelling SL 16 Fr. Voice Prosthesis	05,07 und 09 mm	
IN16xx-SO	CLASSIC Indwelling STERIL 16 Fr. Voice Prosthesis	08-14 mm	
Introduction with gel capsule technique			
IN20xx-NS	CLASSIC Indwelling 20 Fr. Voice Prosthesis	06-20 mm	
AD20xx-NS	ADVANTAGE Indwelling 20 Fr. Voice Prosthesis	06-14 mm	
AS20xx-NS	ADVANTAGE Indwelling 20 Fr. Voice Prosthesis	06-14 mm	
IN20xx-LF	CLASSIC Indwelling LF 20 Fr. Voice Prosthesis	06-14 mm	
IN20xx-LEF	CLASSIC Indwelling LEF 20 Fr. Voice Prosthesis	06-14 mm	
IN20xx-SL	CLASSIC Indwelling SL 20 Fr. Voice Prosthesis	05,07 und 09 mm	
IN20xx-LEIR	CLASSIC Indwelling LEIR 20 Fr. Voice Prosthesis	05 und 07 mm	
IN20xx-LESL	CLASSIC Indwelling LESL 20 Fr. Voice Prosthesis	05 und 07 mm	
IN20xx-IR	CLASSIC Indwelling IR 20 Fr. Voice Prosthesis	06-12 mm	
IN20xx-LETO	TEP Occluder Indw. LETO 20 Fr. Voice Prosthesis	06-14 mm	



Andreas Fahl Medizintechnik-Vertrieb GmbH is the exclusive distributor for Blom-Singer® products in selected European countries. Blom-Singer® is a registered trademark of Freudenberg Medical, LLC.