




# LARYVOX® EXTRA HME SPORT

- ▶ HME in discrete design with speech option via easy finger lock
- ▶ Due to its very low respiratory resistance, the LARYVOX® EXTRA HME SPORT is well suited for use in sporting activities
- ▶ Excellent tactility: a rounded recess in the lid of the HME simplifies the locking mechanism and makes optimal positioning of the finger very easy to feel
- ▶ The lateral openings of the LARYVOX® EXTRA HME SPORT allow for largely unimpaired and easy breathing – even below clothing.
- ▶ A rabbet around the housing allows easy and quick removal of the HME (e.g. when coughing)
- ▶ With 22 mm adapter (KOMBI) to connect to any conventional base plates (e.g. LARYVOX® Tape) or silicone tubes for laryngectomees (e.g. LARYNGOTEC®)
- ▶ In addition to the LARYVOX® EXTRA HME SPORT, the following additional variants are available:

**Normal** - normal breathing resistance at normal physical exertion

**Medium** - medium breathing resistance at medium physical exertion

**HighFlow** - reduced breathing resistance at pronounced physical exertion

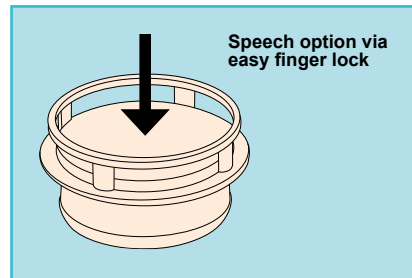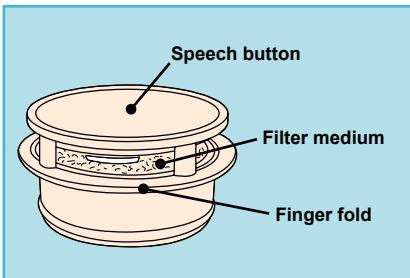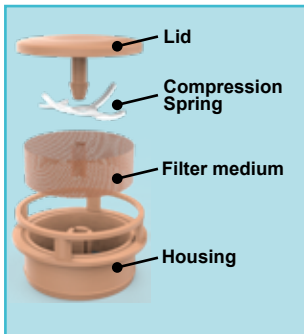

## Order information:

LARYVOX® EXTRA HME SPORT

Pack of 30

REF 49863

LARYVOX® EXTRA HME

Pack of 30

REF 49860

LARYVOX® EXTRA HME MEDIUM

Pack of 30

REF 49862

LARYVOX® EXTRA HME HighFlow

Pack of 30

REF 49861

